



AGELESS
WONDER™



**FACIAL MUSCLE
STIMULATION SYSTEM**

USER'S GUIDE

TROUBLESHOOTING GUIDE

PROBLEM	POSSIBLE CAUSE	SOLUTION
LCD Display does not light up at all	<ul style="list-style-type: none"> - Batteries not installed - Batteries installed incorrectly or battery tray not properly locked into place - Batteries are dead 	<ul style="list-style-type: none"> - Install Batteries - Check polarity "+", "-" - Remove the battery tray and re-insert it until it locks into place. A confirmation beep should sound. - Replace Batteries
LCD Display lights up, but no sensation is felt.	<ul style="list-style-type: none"> - No skin contact - Application wand is not plugged in correctly - Application wand is not attached to headset correctly - Microprocessor freeze - Battery is low 	<ul style="list-style-type: none"> - Check to see if the conductive sponge is wet to the touch - Check that wand plug is connected firmly to the control unit - Check operation with handheld wand first to isolate the problem. If handheld operation works, attach application wand to the headset and- - Check that both wands are properly attached to the headset. - Choose mode and Increase intensity by pushing the + button - Perform a RESET, remove batteries for 1 minute for full RESET - Change Batteries
LCD Display lights up, but nothing occurs.	<ul style="list-style-type: none"> - Intensity set to low level - Conductive sponge is not wet enough - The wands are positioned incorrectly - Conductive sponge has oil deposits. - Battery is low 	<ul style="list-style-type: none"> - Choose mode and Increase intensity by pushing the + button - Put more water on the conductive sponge - Move wand around to locate your trigger points. - Clean or replace conductive sponge. - Change Batteries.
LCD Display lights up, but none of buttons work.	<ul style="list-style-type: none"> - Microprocessor freeze - Battery is low 	<ul style="list-style-type: none"> - Perform a RESET, remove batteries for 1 minute for full RESET - Change Battery

IMPORTANT: POWER RESET

If the unit does not switch on, or if you can not change to different programs, the microprocessor may have "frozen". Please perform a simple **RESET**. Remove batteries and replace them after 1 minute.

THANK YOU

On behalf of the entire team at Thane Direct, we appreciate your recent purchase of AGELESS WONDER™.

Ageless Wonder™ delivers gentle microcurrent impulses to minimize wrinkles and to restore and tone facial muscles for a more youthful appearance.

- Improve facial tone
- Rejuvenate the skin
- Reduce puffiness
- Even out skin tone

Enjoy your AGELESS WONDER™
and, once again,
thank you!

CONTENTS

IMPORTANT!

Please read the complete instruction manual before you use the AGELESS WONDER™. It will give you a better understanding of how the device works and you will be able to enjoy the best results.

Warning	3
Caution	4
Introduction	5
Components	6
Operation	7 - 13
Cleaning and Care	14
Frequently Asked Questions	15
Troubleshooting	16

FREQUENTLY ASKED QUESTIONS


- Q. How does the AGELESS WONDER™ work?
- A. The AGELESS WONDER™ is designed to copy the messages sent from your brain to the facial muscles. These messages cause the muscles to contract and relax.
- Q. Can muscles become over-developed or strained?
- A. No. The principle is the same for general training of any muscle—the more exercise with suitable periods of rest, the stronger and more toned the muscles will become.
- Q. Are there any age restrictions on use of the AGELESS WONDER™?
- A. The AGELESS WONDER™ is for adult use, The elderly should use the device only under supervision and at low intensity. The device should not be used by infants, toddlers or children.
- Q. Will I find it uncomfortable?
- A. As long as the conductive sponge is moist and kept firmly against the skin, you will not feel any unpleasantness. At low intensity, you will feel a slight tingling sensation, but as the intensity gets higher you will only feel pleasant contractions. However; should you feel irritation or uncomfortable sensations continue in a specific area and move on a different area on your face.
- Q. Should the intensity settings be the same for each muscle?
- A. Each muscle will respond differently, so intensities will vary accordingly. Move the application wand around to find the right trigger point that gives you maximum contractions.
- Q. Are there any contra-indications to use of the AGELESS WONDER™?
- A. Please read all the warnings and cautions page before you use.

CLEANING AND CARE


Never immerse the device, the application wand or the headset in water!

Replacing the conductive sponge on the application wand:

- Take the application wand and remove the top plastic cover
- Place the conductive sponges onto the wand and press to replace the plastic top cover
- Use the provided bottle to drip water until the conductive sponge is wet to the touch.
- Clean accessories with a slightly wet sponge with detergent and wipe with a cloth.
- Clean the conductive sponge by soaking it in water with disinfectant diluted in water at 1/10th strength.

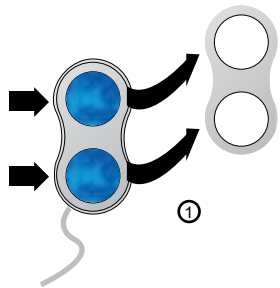

1. Open the plastic top cover of the application wand as shown above using the dent on the side of the application wand.

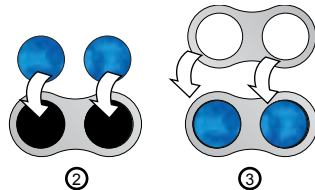

2-3. Replace your old conductive sponges with the new once included in the box and put back the top plastic cover, press to insure it is placed firmly.

WARNING

The AGELESS WONDER™ is extremely safe and simple to use. We assume no responsibility for any harm caused by misuse of the product. Please read the entire instruction manual before using the device.

The AGELESS WONDER™ is intended only for cosmetic use on facial tissue for the stimulation of facial muscles as indicated in the instruction manual for personal beauty purposes. All other uses shall be deemed improper. The long-term effects of electrical stimulation are unknown.

- Use AGELESS WONDER™ only on normal, healthy and clean skin.
- Do not apply stimulation:
 - Over swollen, infected or inflamed areas of skin or skin eruptions.
 - Over or in proximity to cancerous or other skin lesions or open wounds.
 - After experiencing acute trauma or a recent surgery.
 - Through the head.
 - Over your chest as the introduction of electrical current into the heart may cause rhythm disturbances, which could be harmful.
- Never use the AGELESS WONDER™
 - while driving, operating machinery or during activities in which involuntary muscle contractions may endanger the user or others,
 - while sleeping,
 - in the bath or shower.
- Never immerse the AGELESS WONDER™ device or application wand in water or any liquid.
- Keep the AGELESS WONDER™ out of the reach of children.

CAUTION


- Do not use this product if you have a tendency to hemorrhage.
- Do not use this product in areas on your face in which you have injections (such as Botox or facial dermal fillers), or near areas in which you have had surgical implants.
- Do not use over skin areas that lack normal sensation.
- Persons with heart problems should consult their doctor before using product.
- The elderly should use the device only under supervision and at low intensity.
- Immediately stop using this device if you experience any adverse effects such as skin irritation, hypersensitivity, burns or any sort of muscular irritation or discomfort.
- Pregnant women should consult their doctor before using this product.
- Only use the application wand in accordance with the illustrations in this manual.
- Users should start all sessions in a seated position.
- If you have the following conditions, please consult your doctor before using the device:
 - Heart disease or cardiac infirmities, whether or not you have a pacemaker or other implanted device
 - Epilepsy
 - Multiple Sclerosis
 - Phlebitis in its active phase (inflammation of a vein)
 - Recent scar tissue or surgery
 - Varicose veins in later stages

TIMER



Timer

Default time is 5 minutes. Press the **+** / **-** button once to increase or decrease the time increment. Maximum time is 20 minutes.

1. To switch the device off, press and hold the **ON/OFF**  button for 2 seconds.
2. To turn the sound on/off, press and hold the **MODE M** and “-” button for 2 seconds. The display will show the ‘sound off’ icon.

Auto Functions

The device will automatically turn off after 30 seconds if not in use. The **AGELESS WONDER™** device has an auto-timer. It will automatically switch off after the chosen timer setting. If you change the program mode during use, the timer will continue to countdown from the original leftover time and will not restart timer.

What to do if you do not feel any muscle contractions

If no muscle contractions are felt, move the application wand slowly around the area you want to target. Make sure the conductive sponge of the application wand is moist to the touch. Gradually increase the **PULSE INTENSITY** by pushing the **+** button until muscle contraction is felt.


The **AGELESS WONDER™** is a great way to tone, lift and firm your face. It may take a few moments of practice to locate the application wand correctly and apply the ideal pressure that will produce the best results for you.

Tip! *If the conductive sponge of the application wand is not moist enough or the skin contact is not sufficient, you may experience an uncomfortable tingling sensation. The same may happen if the conductive sponge is worn and needs to be replaced.*

INTENSITY



Power with 30 Intensity Levels

Push the **+ / -** button when in the Power Mode  intensity level. Level 1 is the first intensity you can feel. The pre-set value is 0 after choosing your program in the mode setting. After increasing the intensity, you will notice a tingling sensation prior to the muscle contracting. Once the muscle is contracting, you will feel pleasant and relaxed. Adjust the intensity according to your comfort level.

The intensity has 30 levels. With every push of the **+ / -** button, you increase / decrease the intensity by one level.

Press and hold the **+ / -** button for 1.5 seconds to enter fast scroll mode.

Tip! *For best results, we recommend to set the intensity at a low level. To find the right level of intensity, increase the intensity to the maximum tolerable level and reduce it 1-3 levels until you feel comfortable and the unit is still generating strong muscle contractions.*

INTRODUCTION

The AGELESS WONDER™ is a battery-powered portable EMS device intended for the improvement of facial muscle tone and for the tightening and firming of the facial skin. It uses specially designed bi-polar low voltage micro-current impulses to stimulate your facial muscles into contracting.

The face is one of the hardest muscle areas to train and work out. With the AGELESS WONDER™, you can contract the muscle fibers of your face safely and effectively to exercise your facial muscles so they become firmer and give your face a more lifted appearance.

The AGELESS WONDER™ generates small, gentle electrical pulses that are delivered onto your skin through the application wand. These electrical pulses activate your underlying motor nerves, resulting in muscle contractions. The AGELESS WONDER™ sends a signal to your muscle to cause it to contract and then relax. These impulses are designed to mimic the natural nerve signal and to obtain efficient, pleasant, continuous contractions. It stimulates these impossible-to-exercise facial muscles and it is easy to use!

Clinical research trials have been conducted to test the effectiveness of the AGELESS WONDER™ device. The performance of the AGELESS WONDER™ device has been enhanced with 6 pre-programmed modes for ease of use with maximum benefits and results. The programs are specifically designed for 5 different facial areas.

COMPONENTS



- AGELESS WONDER™ control panel
- Headset
- Application wand
- Application wand with connecting wire
- Application wand Handle

- Dripping water Bottle
- X 6 Conductive replacement sponges
- User Guide
- 2 X AAA Batteries

Use the handheld wand for this stimulation mode. Use the wand attached to the cable and plug it into the device. Please note that without the headset, only the application wand with connecting wire will work. We recommend clipping the provided handle into the back of the application wand.

MODE 4 - Forehead

- Select this mode to stimulate the forehead and temporal region.



MODE 5 - Chin, jaw line

- Select this mode to stimulate the chin and lower jaw line.



MODE 6 - Under eye area

- Select this mode to stimulate the area under your eye.



PROGRAM MODES AND APPLICATION AREAS

The AGELESS WONDER™ works on facial muscles that have never been stimulated before. Consistency is the key. Use your AGELESS WONDER™ for 15 minutes a day for 4 weeks. Use the headset for this stimulation mode.

MODE 1 – Cheeks - Lift

- Position the headset horizontally with the application wands below the cheekbone and approximately 1-2 cm in front of the ear.

This program will act to give you a more lifted appearance and also gives your skin a tighter and smoother look and feel.



MODE 2 – Cheeks – Firm and Tone

- Position the headset horizontally with the application wands below the cheekbone and approximately 1-2 cm in front of the ear.

This program will act to give your face a firmer and more toned appearance and also gives your skin a tighter and smoother look and feel.



MODE 3 – Lower Cheek


- Position the headset horizontally with the application wands below the lower cheekbone and approximately 2-3 cm in front of the ear.

Select this mode to stimulate the lower cheek, from the corner of the mouth to the upper jaw.




OPERATION

LCD DISPLAY

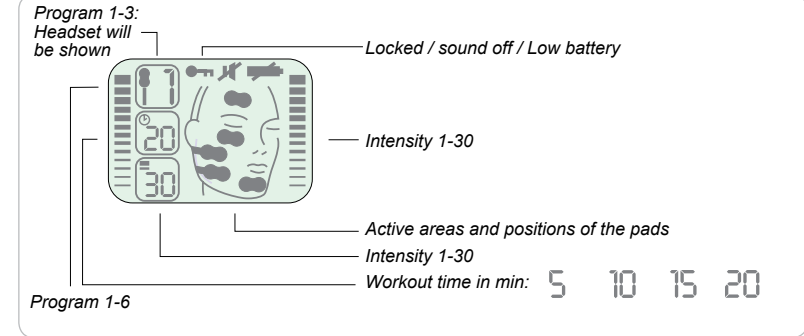


HOW TO OPERATE CONTROLS

Indication:

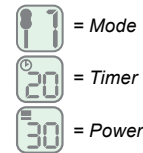

Controls:

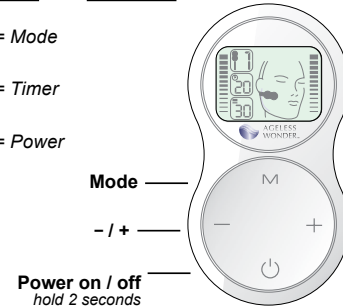

- KEY LOCK TOGGLE ON/OFF:** Press and hold the **MODE M** button together with the “ + ” while the device wand is operating and you are in intensity mode. The LCD will show the key lock icon.
- MUTE TOGGLE ON/OFF:** Press and hold the **MODE M** button together with the “ - ” button for 2 seconds. The LCD will show the mute icon.

Low battery indicator (The AGELESS WONDER™ requires 2 x AAA batteries)

Before using your AGELESS WONDER™.


1. Open the battery compartment by sliding the cover out. Pull out the tray and place 2 AAA batteries in the proper position. Push the tray back in and replace the battery compartment cover.



2. Using the bottle provided, apply water to the conductive sponges of the application wand until they are wet to the touch.



3. The application wands are flexible to follow the contour of your face. Practice with the headset ease of use.


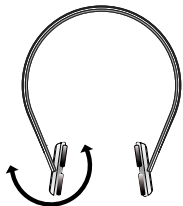


Tip! For program mode 1 and 2, the correct positioning is directly under your cheek bone a towards the back of your face.

The Ageless Wonder™ is used with a headset or with a handle. We suggest using the headset for program modes 1-3 (follow Option A). For program modes 4-6 use the handle and refer to Option B for installing.



Option A: Headset Use (for program modes 1-3):

Clip the application wands to each end of the headset. Once the application wands are securely installed on the headset, plug the wire attached to one of the wands into the unit, turn the unit ON and select your options.

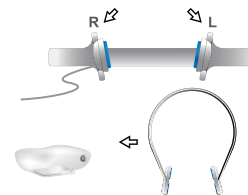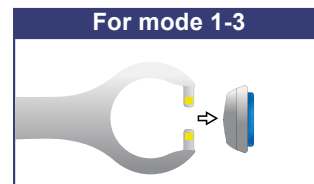

Option B: Handle Use (for program modes 4-6):

Clip the application wand that has wire attached to it into provided handle. Once the application wand is securely installed on the handle, plug the wire attached to the wand into the unit, turn the unit ON and select your options.

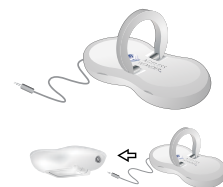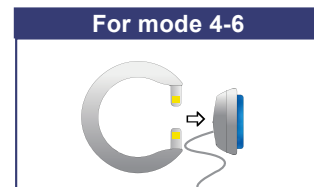



AGELESS
WONDER™



THANE

Distributed by:

Thane International, Inc.
La Quinta, CA 92253, U.S.A.
www.thane.com

Thane Direct Canada Inc.
Mississauga, ON L4W 5M6, Canada
www.thane.ca

Thane Direct UK Ltd.
WC2A 3LH
www.thanedirect.co.uk

Danoz Direct Pty Ltd
Unit D 23-27 Bourke Road
Alexandria, Sydney
NSW 2015, Australia
www.danozdirect.com.au

Copyright © 2011 - 2012 Thane International, Inc.